What Does It Measure?

The Skin-pH-Meter PH 905 is a quick, easy and economical tool to specifically measure the **pH on the skin surface** or the scalp.

The Measuring Principle

The measurement is based on a high quality **combined electrode**, where both H⁺ ion sensitive electrode and additional reference electrode are placed in one housing. It is connected to a probe handle containing the measurement electronics.


Fields of Application

There are many fields of application where changes in skin pH are of interest.

- In cosmetological and pharmaceutical application fields for the development of home and personal care products.
- The role of the pH has gained importance in skin health and is therefore subject to basic research.
- Studies on the changes in the microbiotic environment on skin and scalp.
- For educating on skin health and suitable products in **occupational health**.

Advantages

- The modern, high quality electronics of the probe allow a very quick (1 s) and reliable measurement avoiding occlusion effects.
- The probe head is planar for measuring on the skin surface.
- Continuous measurements over long-term periodes possible, e.g. for in-vitro applications.
- Display of the pH-value with one decimal.
- Regular calibration can easily be done by the user.
- AvailableforC+K MPA-systems, as stand-alone device (MDD) and wireless probe (operation with RR 200 & MPA WLplus software).

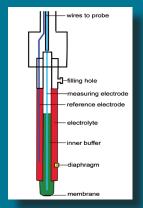









Technical Data

Dimensions: 26.5 cm, Measuring surface: Ø 12 mm flat; Weight: 116 g Measurement range: pH 1,0 to pH 11,0, Measurement uncertainty: \pm pH 0.1 Technical changes may be made without prior notice.

Courage+Khazaka electronic GmbH since 1986 Mathias-Brüggen-Str. 91 · 50829 Köln · GERMANY

phone +49 221 95 64 99 0 · fax +49 221 95 64 99 1 info@courage-khazaka.de · www.courage-khazaka.de

